




BWS-M-G Series

Automatic Bottle Filling Machine

Automatic Bottle Filling Machine



Dedicated to animal drinking water bottle filling, with high efficiency, high degree of automation, simple operation, etc. Can completely replace manual filling.

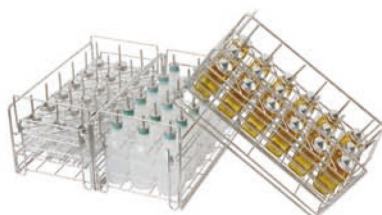

- High Efficiency: Can fill 36 drinking bottles one time (about 1 min).
- Online disinfection: With automatic disinfecting device, which can realize automatic disinfection of piping system, and eliminate pipe pollution.
- Easy operation: To set parameters according to the different volume .
- With automatic quantitative charging device, can add minerals to water.
- Piping design: Fast disassembly type sanitary piping design, convenient for pipeline cleaning and disinfection.
- Filling homogeneity: Each water bottle filling amount uniform, error is not more than 5%.
- Wide scope of application: According to the different volume water bottles, choose a different basket.



Specifications

Code	Name	Size (L x W x H)mm	Weight	Power
BWS-M-G360 BWS-M-G180	Automatic bottle filling machine	770 x 750 x 1530	100	220V 50Hz 200W
Parameters	Filling (bottle/cycle)	Pure Water Pressure /Mpa	Water Inlet	Water Outlet
	36/18	0.2-0.4	3/4"	1/2"

Developed Assembled & Manufactured by

ATNT Laboratories in technical collaboration with

Shinva Medical P.R. China

ISPA engineering Pvt Ltd India

Wagle Industrial Estate, Mumbai

Thane -400604 Maharashtra India

ATNT Laboratories

Unit # 812, Excellencia Lodha Supremus 1,
Rd No.22, Wagle Estate, Mumbai, Thane - 400604, India.

Mob. No. +91 - 9892520959

Tel. +91-22-25830958 / + 91-22-25830959

Email ID : info@atntlabs.com, ashutosh@atntlabs.com

Web : www.atntlabs.com