




## **BSE - CA Series Bedding Disposal Instrument**

## Bedding Collection Instrument



The BSE-CA series bedding collection table is specially used for the animal lab. It is used to prevent the staff from inhaling dust or polluted gas when dumping the cage bedding, and protect the safety of the staff.

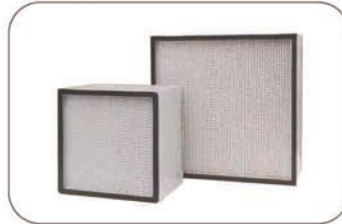

- Operating area for negative pressure, effectively adsorption of dust and odor, protect the safety of personnel.
- Primary effect/high efficiency double filtration, the polluted gas can be directly discharged into the room or connected with pipes to discharge the outdoor air after filtering by HEPA (99.995% high efficiency filter) to improve indoor air quality.
- SUS304 stainless steel dumping table, electrostatic spraying outer cover, easy to clean.
- The dumping port is horn shaped and inclined at a large angle, which is convenient for material collection.
- The dumping port is equipped with collision bar, cage box knock against collision bar, easy to remove silt pad material.
- There are large glass Windows on both sides of the equipment and lights on the top to increase light transmittance.
- Handrails are installed on both sides with casters at the bottom for easy movement.
- The bottom of the equipment is equipped with a large number of movable dirt barrels, which is convenient to collect waste mat bedding.



name	operator	Model	Overall dimensions ( L*W*H )	workspace dimensions ( L*W*H )	The collecting barrel
Bedding collection table	single	BSE-CA-A1000	1100 x 880 x 1950	1000 x 600	2

### Developed Assembled & Manufactured by

ATNT Laboratories in technical collaboration with

Shinva Medical P.R. China

ISPA engineering Pvt Ltd India

Wagle Industrial Estate, Mumbai

Thane -400604 Maharashtra India

### ATNT Laboratories

Unit # 812, Excellencia Lodha Supremus 1,

Rd No.22, Wagle Estate, Mumbai, Thane - 400604, India.

Mob. No. +91 - 9892520959

Tel. +91-22-25830958 / + 91-22-25830959

Email ID : info@atntlabs.com, ashutosh@atntlabs.com

Web : www.atntlabs.com